IMPORTANT SAFETY INSTRUCTIONS

**READ ALL INSTRUCTIONS
BEFORE USING SKN BY
CONAIR™ PORE PURIFIER
ADVANCED MICRO-
DERMABRASION TOOL**

WARNING — To reduce the risk of burns, electrocution, fire, injury or death to persons:

- Plug immediately and discontinue use.
 Not for use by children.
 Do not pull, twist or wrap the cord around the
 appliance.
 Never drop or insert any object into any opening
 in the appliance.
THIS PRODUCT IS FOR HOUSEHOLD USE ONLY.
 Do not leave appliance unattended when plugged
 in.
 Unplug from outlet when not in use and before
 cleaning and servicing.
 Do not unplug by pulling on cord. To unplug,
 grasp the plug, not the cord.
 Do not handle plug or appliance with wet hands.

1. Dry location use only. Do not expose to liquids, vapor, or rain.
2. Ensure jack is plugged into the tool before plugging the USB cable into the outlet.
3. Do not operate appliance with a damaged housing use charging cord that has been damaged.
4. Always unplug the unit before cleaning or when not in use.
5. Be sure voltage is correct before plugging in. The rated input is 5V DC 1.0A.

SAVE THESE INSTRUCTIONS

The SKN by conair™ pore purifier advanced microdermabrasion tool is designed to gently clean your pores and exfoliate your skin without irritation, for healthy-looking results.

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TWO SPEED OPTIONS

The two speeds available on your SKN by conair™ pore purifier advanced microdermabrasion tool allow you to customize your skincare routine to suit your individual needs.

The warm plate on the bottom of the tool helps to open pores, leaving skin feeling bright and smooth.

SMALL PORE EXTRACTOR: To remove blackheads & whiteheads.

2 CRYSTAL DIAMOND HEAD TIPS: For deep exfoliation.

The SKN by conair™ pore purifier advanced microdermabrasion tool is cordless and requires a full charge before use. After tool is fully charged, it can provide up to 100 minutes of use without recharging.

Do not throw device into fire, apply heat to it, or charge, use in a high-temperature environment.

1. Before using the SKN by conair™ pore purifier advanced microdermabrasion tool for the first time, it must be fully charged for 2.5–3 hours. Plug the DC jack into the side of the unit, then plug power cable and adapter into wall outlet. The white LED will begin to flash, indicating that the tool is charging.

2. The charging indicator will FLASH WHITE when charging. Once fully charged, the light will illuminate SOLID WHITE. Always check to be sure that the charging light illuminates when you plug the charging cable into the outlet.

3. The white light will start flashing for 5 seconds and then shut off when the unit reaches low battery capacity during operation.

4. The adapter may become warm while charging — this is normal. Once fully charged, the device can provide 100 minutes of use.

5. In order to maintain the optimum capacity of the rechargeable battery, the SKN by conair™ pore purifier advanced microdermabrasion tool should be given a full 2.5–3 hour charge every 3 months.

6. Remove charging cable from unit and unplug adapter from outlet when not in use. Do not unplug the adapter from the wall when it is still connected to the unit. This will cause the battery to discharge.

Warm Plate and Attachments

- NOTE: Do not do more than two passes over the same area during treatment.**

**STARTING THE SKN BY CONAIR™ PORE
PURIFIER ADVANCED MICRODERMABRASION
TOOL AND SELECTING THE SPEED/SUCTION
LEVEL**

Press and hold the Power Button to turn on the device. The bottom light will illuminate and the WARM PLATE will heat up instantly.

Press the Power Button one more time to start low
 suction setting.

For a deeper exfoliation or stronger suction, press the Power Button again.


To stop the suction at any time and switch back to the WARM PLATE only, press the Warm Plate Button

At any time, the SKN by conair™ pore purifier
advanced microdermabrasion tool can be stopped by
pressing and holding the power switch button
for more than 1 second.

Press and hold the Power Button for 1 second to activate the WARM PLATE. After 5 minutes in this setting, the device will go into standby mode.

At any time, the WARM PLATE can be stopped by pressing and holding the Power Button for more than 1 second.

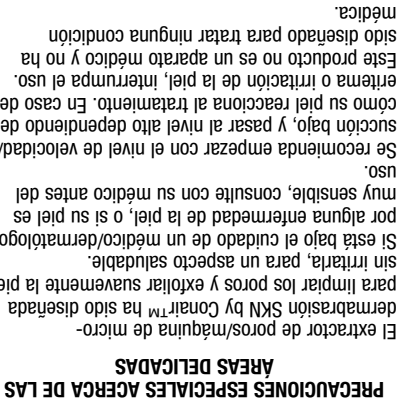To attach the different attachments, hold the unit and push the attachment onto the post until it snaps in place. To remove attachment head, simply pull attachment off of the unit.

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